Withings Scan Monitor 2.0 with Withings Body Scan

Withings Scan Monitor 2.0 Product Guide

Disclaimer: Information in this guide may change without notice.

R_x ONLY: Federal law (USA) restricts this device to sale by or on the order of a licensed practitioner.

WITHINGS

TABLE OF CONTENTS

Sommaire

ENGLISH	03
FRANÇAIS	17
SPANISH	31

PRODUCT GUIDE

This guide explains how to use your Withings Scan Monitor 2.0 functionalities on your Withings Body Scan.

It applies to users in the United States of America.

Intended Use

The Withings Scan Monitor 2.0 is intended to record, store and transfer lead II and lead-III of a two-channel electrocardiogram (ECG). In addition, it calculates and displays leads I, aVR, aVF, aVL.

The Withings Scan Monitor 2.0 also displays ECG rhythms and detects the presence of atrial fibrillation (when the monitor is prescribed or used under care of a physician). The Withings Scan Monitor 2.0 is intended for use by healthcare professionals, patients with known or suspected heart conditions and health-conscious individuals.

The Withings Scan Monitor 2.0 is intended for adult patients in hospitals, clinics, long-term care facilities, and homes.

Contraindications

- **DO NOT USE the Withings Scan Monitor 2.0 with a cardiac pacemaker, ICD, or other implanted electronic devices.**

- **DO NOT USE if diagnosed with an arrhythmia other than Atrial Fibrillation or diagnosed with bundle branch block.**

- Withings Scan Monitor 2.0 is not intended to continuously monitor vital signs in critical conditions or where the nature of the variations is such that it could result in immediate danger to the patient. Withings Scan Monitor 2.0 does not provide alarms.

- Withings Scan Monitor 2.0 is not intended to provide or to replace a diagnosis by a physician or qualified health care professionals. Vital signs measurements, such as those taken with this device, cannot identify all diseases. **If you think you may be having a heart attack (myocardial infarction) or are facing a medical emergency, call the emergency services.**

Warnings & safety notice

- Regardless of the measurements taken using this device, you should consult your doctor immediately if you experience symptoms that could indicate you are experiencing a sudden and/or severe change in health.
- The Scan Monitor Withings 2.0 cannot detect heart attacks or ischemic heart conditions. If you ever experience chest pain, pressure, tightness or what you think is a heart attack, call emergency services.
- **DO NOT** self-diagnose or self-medicate on the basis of this device without consulting your doctor. In particular, do not start taking any new medication or change the type and/or dosage of any existing medication without prior approval.
- Withings Scan Monitor 2.0 cannot detect all instances of Atrial Fibrillation. You should contact your physician if you experience any changes to your health.
- Withings Scan Monitor 2.0 cannot detect arrhythmias other than Atrial Fibrillation. You should notify your physician if you detect possible changes in your health.
- **DO NOT** use to diagnose heart-related conditions.
- **DO NOT** use the Withings Scan Monitor 2.0 if it is damaged. Use of a damaged device could cause patient injury or equipment failure.
- **DO NOT** take recordings when Withings Scan Monitor 2.0 is in close vicinity to strong electromagnetic fields (e.g., electromagnetic anti-theft systems, metal detectors).
- **DO NOT** take recordings when Withings Scan Monitor 2.0 is outside of the operational temperature range (5°C - 35°C) indicated in the Withings Body Scan user manual and humidity range of 20% to 80% relative humidity.
- **DO NOT** take a recording during physical activity.
- **DO NOT** try to repair or modify this equipment yourself. No service or repairs should be performed on the Withings Body Scan other than the maintenance listed in the device's instructions for use.
- The power cord of the charger and the handle cable may cause a strangulation effect due to excessive length. Keep it away from children and pets.
- **DO NOT** use the ECG for morphological analysis.
- For users with Body Mass Index (BMI) ≥ 30 , a higher number of "Noise - Inconclusive" results is to be expected. Repeating the measurement by following the guidelines stated in the section "Recording an ECG" might provide a conclusive classification of your ECG measurement.

Important notice

Before using Withings Scan Monitor 2.0, review the information in this guide. You may also find this guide online at: <https://www.withings.com/guides>.

Withings Scan Monitor 2.0 is intended for use in adults (22 years or older in the U.S.)

Instructions for Use of the Withings Body Scan are described in Withings Body Scan's Product Guide provided with the product.

Retain this documentation for future reference.

Please contact Withings when in need of assistance, setting up, using or maintaining the device or to report unexpected operations or events.

Any serious incident that has occurred in relation to the Withings Scan Monitor 2.0 should be reported to Withings and regulatory authorities in your country of residence.

Please note that the device has not been clinically validated for heart rates higher than 140 bpm.

Using Withings Scan Monitor 2.0 with Withings Body Scan

On-boarding for new Withings' users:

Before being able to use the Withings Scan Monitor 2.0 with the Withings Body Scan, you will need to download the Withings companion app and create an account.

- Only use official app stores to download the Withings app. To ensure that the app is the official Withings space, use the following link: <https://go.withings.com>
- The Withings app is compatible with smartphones or tablets with iOS 14 or later, or with Android 10 or later;
- Use a trusted Wi-Fi network with your Withings app. Do not use a public Wi-Fi network you don't know;
- The mobile application is not intended to be used on a computer. No anti-virus software is needed.

You are now ready to set up the scale Withings Body Scan with the Withings app on your iOS or Android device. Please note that the Withings Scan Monitor 2.0 is only compatible with Withings Body Scan.

Setup:

- Once you have downloaded the Withings app from an official store (App store and Google Play store), open the app — see previous page for more information.
- In the Devices tab, select "Install a device," then select "Scales."
- Select the product "Body Scan."
- Follow the on-screen instructions. You will be prompted to pair your device via Bluetooth.
- After pairing, a tutorial section will be available. Follow the ECG tutorial in order to activate Withings Scan Monitor 2.0. You may find more information regarding this activation page 9 "How to enable ECG Measurement?" section).
- You may exit onboarding at any time by tapping Cancel.

Charging:

- Withings Scan Monitor 2.0 is not intended to be used while charging.
- When in charge, Withings Scan Monitor 2.0 does not perform any measurement.

What is a 6-lead ECG?

An ECG, or electrocardiogram, is the graphical representation of the electrical activity of the heart. With each heartbeat, an electrical wave travels through your heart. This wave causes your heart to contract and pump blood.

In a doctor's office, a standard 12-lead ECG is usually taken by placing electrodes on your arms, legs, and torso. A lead measures the difference of electric potential between a pair of electrodes. Each lead measures this from a different angle.

The Withings Scan Monitor 2.0 uses electrodes from the Withings Body Scan to record a 6-lead ECG. Electrodes are in contact with the hands and feet.

A 6-lead ECG is able to provide information about heart rate and certain heart rhythms: Sinus rhythm and Atrial Fibrillation (AFib). 6-lead ECGs are sometimes prescribed by doctors for people to wear at home or within the hospital so that the doctor can get a better look at the underlying rate and rhythm of the heart.

Recording an ECG


How to take an ECG recording:

1 - Make sure the Withings Body Scan is on a stable surface, and that you activated and enabled the ECG measurement in the Withings companion app (see page 9 "How to enable ECG Measurement?" section).


2 - Grab the handle. Hold the handle in both hands, taking care that the part of the handle marked with an "L" is in your left hand and that the part of the handle marked with an "R" is in your right hand. Make sure both hands touch the metallic electrodes on the back of the handle. The handle should be held with the cable positioned vertically, as shown in the figure below.

3 - Step on the scale with bare feet. Make sure your feet are centered on the scale's stripes (electrodes).

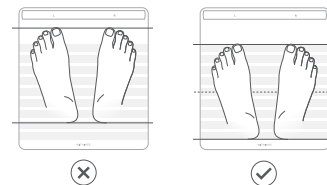4 - Keep your arms along the side of your body and relax. You do not need to grip the handle very tightly.

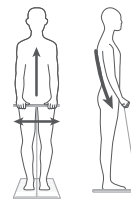5 - A series of measurements will start depending on the configuration you chose in the application (see page 04 "Using Withings Scan Monitor 2.0 with Withings Body Scan"). You will know that the ECG measurement has begun when the ECG recording screen is displayed on the scale LCD screen and after the first vibration. Do not speak or move during the ECG recording.



- It takes 30 seconds to complete a recording.
- Wait for the measurement to end. A countdown indicates the remaining time.
- When the measurement is complete you will feel a double vibration.

ECG Outputs

The Withings Scan Monitor 2.0 utilizes Lead I for classifying an ECG output.

ECG

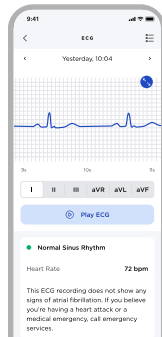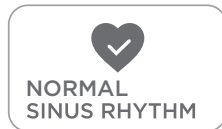Normal Sinus Rhythm

Normal Sinus Rhythm means your heart rate is between 50 and 100 beats per minute (bpm) and is beating in a regular pattern. The recording does not show any signs of atrial fibrillation. If you believe you're having a heart attack or a medical emergency, call emergency services.

ECG



Atrial Fibrillation

This recording shows signs of atrial fibrillation, and your heart rate is beating at least at 50 beats per minute (bpm). Get in touch with your doctor. Atrial Fibrillation occurs when the two upper chambers of the heart move randomly instead of pumping regularly. This does not allow for complete emptying of the chambers and thus, blood may become stagnant and create blood clots.

ECG Outputs

ECG

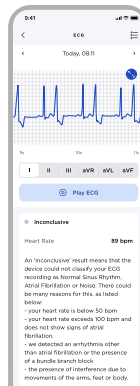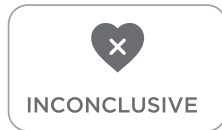
Inconclusive

An 'Inconclusive' result means that the device could not classify your ECG recording as Normal Sinus Rhythm, Atrial Fibrillation or Noise.

There could be many reasons for this, as listed below:

- your heart rate is below 50 bpm
- your heart rate exceeds 100 bpm and does not show signs of atrial fibrillation.
- we detected an arrhythmia other than atrial fibrillation or the presence of a bundle branch block
- the presence of interference due to movements of the arms, feet or body.

ECG

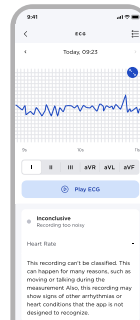Noise - Inconclusive

This recording can't be classified. This can happen for many reasons, such as moving or talking during the measurement. Also, this recording may show signs of other arrhythmias or heart conditions that the app is not designed to recognize.

How to enable ECG measurement?

- First, make sure your Withings Body Scan is configured properly. Then please activate Withings Scan Monitor 2.0 by following the ECG tutorials in the Withings companion app (see below).

- You can activate Withings Scan Monitor 2.0 by going to “Devices” > “Body Scan” > “Tutorials”, then selecting the user you are willing to activate the ECG for, and then following the information on the companion application.

- The first ECG recording will need to be reviewed by a healthcare professional. It should take less than 24 hours. Once completed, the feature will be fully unlocked, and the classification of the Atrial Fibrillation detection algorithm will be displayed both on the device and in the app.

Please take into consideration that the scale will need to be synchronized once before the change is applied to your device. Performing a measurement will apply the changes to your device. If you don't perform a measurement, the changes will be applied to your device during a routine synchronization overnight.

Sharing ECG results

How to share your ECG with your doctor:

Once the ECG measurement feature is activated (see page 04 “Using Withings Scan Monitor 2.0 with Withings Body Scan”), the results of the ECG can be shared with a physician as a PDF, generated by the Withings app. The pdf includes the following information:

- The ECG strip for the 6 Leads
- The average heart rate
- The classification of the ECG according to the detection algorithm.

Troubleshooting – ECG

If you experience difficulties with operating your Withings Scan Monitor 2.0, refer to this troubleshooting guide. If you cannot fix the problem using the troubleshooting instructions, please contact Withings.

- I cannot get the Withings Scan Monitor 2.0 to take an ECG reading.

- Solutions:
- Ensure that you have completed all of the setup steps in the Withings app on your smartphone.
 - Ensure that your legs, feet, arms, hands and fingers remain still during recordings and that your feet and fingers are in contact with the electrodes.
 - Ensure that no rings or other jewelry are touching the electrodes.
 - Make sure the glass plate and handle of your Withings Body Scan are clean and dry. Make sure that your hands and feet are dry. Water and sweat can be the cause of a poor quality ECG.

- I have an inconclusive measurement. It looks like the ECG recording has a lot of artifacts, noise, or interference.

- Solutions:
- Ensure that your legs, feet, arms, hands and fingers remain still during recordings and that your feet and fingertips are in contact with the electrodes.
 - Try to relax and not to move.
 - Make sure the glass plate and handle of your wrist and your Withings Body Scan are clean and dry. Make sure that your hands and feet are dry. Water and sweat can be the cause of a poor quality ECG.
 - Move away from any electronics that are plugged into an outlet to avoid electromagnetic interference.
 - Remove any ring you might be wearing.

- The ECG waveforms appear upside down.

- Solutions:
- The handle may be held with the left electrode in the right hand and the right electrode in the left hand

Technical specifications

Product Name: Withings Scan Monitor 2.0

Model: WBSSM2

ECG Sensors: 2 stainless steel electrodes (handle) and 1 ITO electrode (glass plate)

Battery Operated: 1 year typical use on a single charge

Power Source: Power internally by a 3.6 Vdc Lithium ion battery (use the USB-C charging cable included) and a DC 5V power adapter with a minimum power of 3W and a maximum power of 5W.

Heart Rate Measurement Range: 30 bpm to 220 bpm

Display: Color 3.2inch LCD screen

Operation time: 30 seconds

Transport and storage conditions: Max altitude: 2000m

🌡️ -10°C (14°F) — 50°C (122°F)

💧 10% — 90%

📶 86kPa — 106kPa

Operation conditions:

🌡️ 5°C (41°F) — 35°C (95°F)

💧 20% — 80%

📶 86kPa — 106kPa

IP Protection level : IP22

Mode of operation : Continuous Operation

Minimum product life: 3 years

Please be aware that the stated 3-year product life of Withings scales represents a minimum lifespan under standard tested conditions. Withings products are designed to last: our scales are man-

ufactured with care, and are provided with continual over-the-air firmware updates to avoid obsolescence. This is why with proper care and usage, our health products are likely to serve you reliably for much longer than the minimum product life. In fact, we're proud to share that 50% of our users are still actively using their Withings scales after 7 years.

Wireless information

Wireless Specifications:

Mode	Frequency Band (MHz)	Maximum Output Conducted Power (dBm)
BT LE	2402-2480	5 dBm
WLAN	2412-2462	16 dBm

The wireless communication of the device is supported by BLE and Wi-Fi communication. This communication is established between the device and the Withings application. The communication between the device and the Withings application is encrypted through an exchange of a paired key.

The communication latency between the device and the Withings application takes less than 10 seconds when the device and the smartphone are less than 16 feet (5 meters) apart.

The communication between the device and the Withings application is not modified with sources of interference located within 5 meters. Wireless coexistence has been tested in accordance with the following standards:

- ANSI C63.27:2017 and,
- AAMI TIR69:2017

Electromagnetic disturbances have been tested in accordance with the standard IEC 60601-1-2:2021.

Cleaning, Disinfection, Maintenance and Storage

Cleaning, Disinfection and maintenance of Withings Body Scan

- It is recommended to clean the glass plate and the handle electrodes of the scale before each measurement.
- Use only a cloth dampened with warm water to clean your scale. The use of any detergents or cleaners may cause irreversible damage to the conductive coating on top of the glass plate.
- [FOR HOSPITAL USE ONLY] For disinfection you may use Liquinox (Alconox). Please follow the manufacturer's instructions for use.
- Do not use the device while doing cleaning and maintenance operations.
- To recharge the battery, follow instructions in this guide.

Storage:

- After a measurement using the handle, lower slowly so it rolls back into the scale without any sudden drop.
- Store the device in a clean and safe location between use.
- If storage conditions are different from the operating conditions indicated in the technical specifications, wait 30 minutes before taking a measurement.

Safety and performance

- Withings Scan Monitor 2.0 software's ability to accurately classify an ECG recording into Atrial Fibrillation and Normal Sinus Rhythm categories was tested according to the IEC 60601-2-47* standard and a clinical study conducted by Withings.
- On five public databases, the Withings Scan Monitor 2.0 demonstrated 99.06% sensitivity in classifying Atrial Fibrillation and 99.40% specificity in classifying normal sinus rhythm in classifiable recordings. In the beat-to-beat detection of QRS complexes, the Withings Scan Monitor 2.0 reached a F1-score of at least 99.19% on all the datasets, with the exception of NSTDB where the F1-score was 90.65% because of digitally added noise.
- In the clinical study conducted by Withings, rhythm classification by Withings Scan Monitor 2.0 was compared to ECG recordings reviewed by cardiologists in a validation study with 321 patients. 23.3% of recordings were inconclusive. On conclusive recordings, the sensitivity in classifying Atrial Fibrillation was 99% (lower bound of the 95% confidence interval: 93%) and the specificity in detecting normal sinus rhythm was 99% (lower bound of the 95% confidence interval: 97%).

*IEC 60601-2-47:2012: Requirements for the Basic Safety and Essential Performance of Ambulatory Electrocardiographic Systems.

Safety and performance

- The table below provides the device performance (confusion matrix) for all 4 outputs of Withings Scan Monitor 2.0 when compared to a 12-lead ECG reference standard:

Withings Clinical Study (321 participants)		12-lead Reference ECG signals reviewed by cardiologists				
		Normal Sinus Rhythm	Atrial Fibrillation	Inconclusive	Noise	TOTAL
6-lead Withings Scan Monitor 2.0	Normal Sinus Rhythm	170	1	9	5	185
	Atrial Fibrillation	1	74	0	2	77
	Inconclusive	4	8	21	1	34
	Noise - Inconclusive	3	18	2	2	25
	TOTAL	178	101	32	10	321

- The table below provides the performance for the Withings Scan Monitor 2.0 compared to a FDA cleared device, also known as a predicate device (Withings Scan Monitor):

Atrial Fibrillation (AF) /Normal Sinus Rhythm (NSR)	Withings Scan Monitor 2.0	Withings Scan Monitor
Sensitivity [95% confidence interval lower bound]	0.73 [0.639]	0.77 [0.675]
Specificity [95% confidence interval lower bound]	0.955 [0.914]	0.965 [0.912]
Sensitivity AF/NSR only [95% confidence interval lower bound] ("Inconclusive" and "Noise" excluded.)	0.989 [0.928]	0.963 [0.894]
Specificity AF/NSR only [95% confidence interval lower bound] ("Inconclusive" and "Noise" excluded.)	0.994 [0.968]	1.000 [0.967]

The sensitivity of Withings Scan Monitor was estimated to be 77% (lower bound 95%CI=67.5%) whereas the sensitivity of the Withings Scan Monitor 2.0 was 73.3% (lower bound 95%CI=63.9%). Although there are differences in sensitivities between the two devices, both devices have equivalent performance because both devices use the same algorithm and meet the confidence intervals criteria and the difference in terms of sensitivity is not statistically significant.

- These results reflect use in a controlled environment. Real-world use of Withings Scan Monitor 2.0 may result in a greater number of strips being deemed inconclusive.

Safety and performance

RF Statement

Medical electrical equipment needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided in the accompanying documents.
Portable and mobile RF communications equipment can affect medical electrical equipment.

Guidance and manufacturer's declaration-electromagnetic emissions

The device is intended for use in the electromagnetic environment specified below. The customer or the user of the The device should ensure that it is used in such an environment.

Emissions Test	Compliance	Electromagnetic environment – guidance
CE emissions CISPR11	Group 1, Class B	The device uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RE emissions CISPR11	Group 1, Class B	
Harmonic emissions IEC 61000-3-2	Not applicable	The device is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Voltage fluctuations/Flicker emissions IEC 61000-3-3	Not applicable	

Safety and performance

Declaration: electromagnetic emissions and immunity for equipment and systems that are not life-supporting and are specified for use only in a shielded location.

Declaration – electromagnetic immunity

Withings Body Scan is intended for use in the electromagnetic environment specified below. The customer or the user of the Withings Body Scan should ensure that it is used in such an environment.

Immunity test	IEC 60601 Test Level	Electromagnetic environment – guidance
Electrostatic discharge (ESD) IEC 61000-4-2	Contact: ±8 kV Air: ±2 kV, ±4 kV, ±8 kV, ±15 kV	Floors should be wood, concrete, or ceramic tile. If floors are covered with synthetic material the relative humidity should be at least 30%.
Radiated RF EM fields IEC 61000-4-3	10 V/m 80 wMHz to 2.7 GHz+ Proximity fields from RF wireless communications equipment	Portable and mobile RF communications equipment should be used no closer to any part of the equipment or system including cables than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Interference may occur in the vicinity of equipment marked with the following symbol:  .
Electrical Fast Transient/Burst IEC 61000-4-4:2012 ED 3.0	N/A	-
Surge IEC 61000-4-5	N/A	-
Conducted RF IEC 61000-4-6	3 Vrms 24 MHz to 80MHz and 6 Vrms in ISM and amateur radio bands between 24 Mhz and 80 MHz	-
Power frequency (50/ 60Hz) magnetic field IEC 61000-4-8	30 A/m 50 Hz and 60 Hz	The power frequency magnetic field should be at levels characteristic of a typical location in a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	N/A	The main power quality should be of the kind used in a typical commercial or hospital environment. If the user of the equipment or system requires continued operation during power main interruptions, it is recommended that the equipment or system be powered from an uninterruptible power supply or a battery.
Common RF emitters (RFID, Wireless Power Transfer, 5G Cellular)	As defined per AIM Standard 7351731, and IEC 61000-4-39	-

Federal Communication Commission interference statement

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. FCC Caution: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment. This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

IMPORTANT NOTE:

Radiation Exposure Statement:

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with a minimum distance of 20 cm between the radiator & your body.

This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter. Country Code selection feature to be disabled for products marketed to the US/CANADA

Disposal

Actuation of European directive 2012/19/EU, for reduction in use of dangerous substances in the electric and electronic device and for garbage disposal. The symbol applied on the device or its packaging means that at the end of its useful life, the product not be disposed of with domestic waste. At the end of the device's useful life, the user must deliver it to a collection center for electric and electronic garbage, or return it to the retailer when purchasing a new device. Disposing of the product separately prevents possible negative environmental and health consequences deriving from inadequate disposal. It also allows the recovery of materials it is composed of to save energy and resources and avoid negative effects to the environment and health. In case of abusive disposal of devices by the user, will be applied administrative endorsements in compliance with current standards. The device and its parts must be disposed of as appropriate, in accordance with national or regional regulations.

Warranty

Withings™ Two (2) Years Limited Warranty - Withings Body Scan

Withings, 2 rue Maurice Hartmann, 92130 Issy-les-Moulineaux, ("Withings") warrants the Withings-branded hardware product ("Withings Product") against defects in materials and workmanship when used normally in accordance with Withings' published guidelines for a period of TWO (2) YEARS from the date of original retail purchase by the end-user purchaser ("Warranty Period"). Withings' published guidelines include but are not limited to information contained in technical specifications, safety instructions or the quick start guide. Withings does not warrant that the operation of the Withings Product will be uninterrupted or error-free. Withings is not responsible for damage arising from failure to follow instructions relating to the Withings Product's use.

Need help?

→ withings.com/support

Equipment Symbol Description

Do not dispose of this product as unsorted municipal waste, take it to electronic recycling.

Read the Instruction Manual before use

Direct Current

Medical device

Medical device manufacturer

Temperature range

Upper and lower limits of pressure

Upper and lower limits of relative humidity

IP22 Ingress of water or particulate matter

Type BF-Applied Part (handle, glass plate)

Do not use if package is damaged and consult instructions for use

Caution

Complies with FCC regulations

Unique device identifier

Serial number

Security

Withings recommends that you add a passcode (personal identification number [PIN]), Face ID or Touch ID (fingerprint) to your phone to add a layer of security. It is important to secure your phone since you will be storing personal health information. Users will also receive additional update notifications on the device via the app, and updates are delivered wirelessly, encouraging rapid adoption of the latest security fixes.

Do not install the device on a smartphone that you do not own.

Do not use a public Wi-Fi network you don't know. Use a trusted Wi-Fi network with your device.

The Withings application is not intended to be used on a computer. No indication of anti-virus software is needed. Only use official app stores to download the Withings application. In case of doubt, use the link go.withings.com.

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Withings Scan Monitor 2.0 | Proactive Health Tracking © 2023 Withings. All rights reserved. Nov. 2023